

# BruxZir<sup>®</sup> *esthetic* NOW

POSTERIOR BLOCK

*"Beauty in a Single Visit"*

NEW!

## The Fully Sintered, Highly Esthetic Zirconia for In-Office Use

- Recommended for posterior crowns demanding high strength and excellent esthetics
- Mills in under 40 minutes\*
- Available exclusively for Glidewell milling systems
- Optimized combination of strength and beauty



PATIENT-PLEASING  
RESULTS



NO OVEN  
NEEDED



POSTERIOR  
CROWNS



HIGH  
STRENGTH

For more information  
[glidewell.com/milling-blocks](http://glidewell.com/milling-blocks) | 844-949-7184

**glidewell.io**<sup>™</sup>  
IN-OFFICE SOLUTION

## BENEFITS

- Lifelike translucency engineered for the posterior
- Deliver restorations right from the mill with no post-processing required
- Easy delivery — conventional cementing with no need to bond
- Preserve your patients' existing dentition with minimal preparation requirements

## SPECIFICATIONS

### MATERIAL

- Pre-shaded, fully sintered zirconia
- Flexural strength: 870 MPa

### IDEAL FOR

- Single-unit posterior restorations, including bicuspid
- Cases with limited interocclusal space
- Restorations with < 0.7 mm margin thickness

### SHADES

- 16 VITA® Classical Shades: A1, A2, A3, A3.5, A4, B1, B2, B3, B4, C1, C2, C3, C4, D2, D3, D4
- 2 Bleach Shades: BL1, BL White

## READY TO DELIVER IN UNDER 45 MINUTES



MILL

+




FINISH/POLISH

+




NO OVEN  
REQUIRED

=

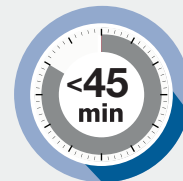

TOTAL TIME

## EASY PREPARATION

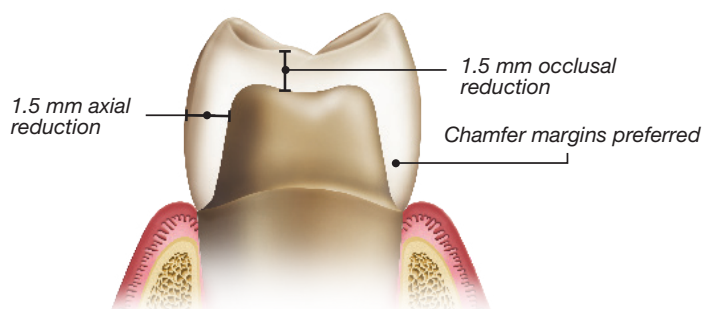

- 0.7 mm incisal/occlusal reduction required (1.5 mm recommended)
- Chamfer or shoulder margins preferred; feather-edge OK
- All axial walls must be convergent

For more information  
[glidewell.com/milling-blocks](http://glidewell.com/milling-blocks) | 844-949-7184

\*Milling times vary based on prescription and mill used. Average milling time shown reflects use of the fastmill.io™ In-Office Mill. VITA is a registered trademark of VITA Zahnfabrik.

 **Glidewell**  
*for the sake of smiles*

PK-1524218-070925